

Life Enhancing. Patient Focused.



Targeted Solutions for: General Imaging, Women's Health, Emergency Medicine, Men's Health, Cardiology, Vascular & Neurosonography

MEDISON

Imaging Accessories - Vol. 2

CIVCO's Needle Guidance Systems utilize a **TWO-PART SYSTEM** consisting of a reusable biopsy bracket and a disposable snap-on needle guide. The disposable needle guides snap onto the biopsy bracket, providing a fixed path for visualizing ultrasound-guided instrument placement procedures. Applications include tissue biopsy, fluid aspirations and catheter placement. **SINGLE-ANGLE** brackets offer imaging professionals a fixed angle needle path during puncture procedures. **MULTI-ANGLE** brackets offer professionals different angles for needle placement. Once the desired angle has been selected, a stainless steel pin locks the angle securely into position. The pin is removable from the bracket for easy reprocessing between patients. Brackets should be cleaned and disinfected according to your CIVCO User's Guide.

STEP 1 **ULTRA-PRO II™ REUSABLE STARTER KITS**

Each Ultra-Pro II starter kit contains 5 sterile procedure kits including: needle guide, telescopically-folded CIV-Flex™ transducer cover, gel packet & colored elastic bands. U.S. Patent No. 5,941,889.

C2-6IC

10cm Depth



Ultra-Pro II™ style bracket
610-888 (Bracket only)
Replacement Kit 610-608

L5-9EC

2.5cm Depth



Ultra-Pro II™ style bracket
612-085
Replacement Kit 610-608

L6-12IS • L8-15IS

1.5, 2.5, 4cm Depth



Ultra-Pro II™ style bracket
678-002
Replacement Kit 610-608
Lock Pin Kit 610-1007


VN6-12 • VNW6-12

2cm Depth



Ultra-Pro II™ style bracket
742-327
Replacement Kit 610-608

Multi-Angle Lock Pin

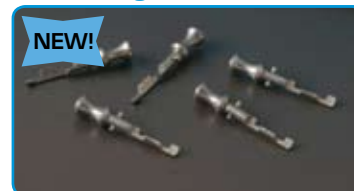

Non-sterile
multi-angle lock pin
610-1007

STEP 2 **ULTRA-PRO II™ DISPOSABLE REPLACEMENT KITS**



ENHANCED PRODUCT! The Ultra-Pro II needle guide now incorporates more flexibility and advanced features to increase productivity during ultrasound-guided procedures. The needle guide directs instruments according to on-screen software guidelines. The Ultra-Pro II features a larger tab for improved quick-release function, allowing easy detachment of the needle from the transducer. Easy-to-read gauge sizes on the inserts makes it simple to identify and alter gauge sizes in a darkened ultrasound suite. The guide now features a larger funnel for instrument insertion and will accept a full range of sizes, including: 8.5FR, 14-23GA (19GA not available). U.S. Patent Nos. 5,941,899.

REF. #	Product Description	Qty.
610-699	Sterile Ultra-Pro II	24
610-579	Sterile Ultra-Pro II with 14 x 61cm (5.5" x 24") CIV-Flex cover (3D)	24
610-608	Sterile Ultra-Pro II with 14 x 91.5cm (5.5" x 36") telescopically-folded CIV-Flex cover (3D)	24
610-609	Sterile Ultra-Pro II with 17.8 x 91.5cm (7" x 36") telescopically-folded CIV-Flex cover (3D)	24
610-904	Sterile Ultra-Pro II with 14 x 147cm (5.5" x 58") telescopically-folded CIV-Flex cover (3D)	24
610-917	Sterile Ultra-Pro II with 17.8 x 147cm (7" x 58") telescopically-folded CIV-Flex cover (3D)	24

STEP 1 MAGGI II Plus™ Reusable Starter Kit

C2-5EL • C5-2EL



8cm Depth

Ultra-Pro II™ style bracket
610-897
Replacement Kit 610-239

STEP 2 MAGGI II Plus™ Replacement Kits



When attached to a reusable bracket, the MAGGI II Plus directs needles, catheters and other instruments in accordance with the ultrasound system's on-screen software guidelines. The MAGGI II Plus accepts instruments 14, 16, 18, 20, and 22 gauge in one convenient insert. The MAGGI II Plus retrofits all existing MAGGI style brackets.

REF. #	Product Description	Qty.
610-237	Sterile MAGGI II Plus	24
610-239	Sterile MAGGI II Plus with 14 x 61 cm (5.5" x 24") flat-folded CIV-Flex cover (3D)	24
610-240	Sterile MAGGI II Plus with 17.8 x 61 cm (7" x 24") flat-folded CIV-Flex cover	24

STEP 1 Multi-Pro 2000™ Reusable Starter Kit

L4-7EL



2.5cm Depth

Multi-Pro II™ style bracket
613-203
Replacement Kit 610-608

STEP 2 Multi-Pro 2000™ Replacement Kits



The Multi-Pro 2000 provides a safe, efficient means of performing ultrasound-guided biopsies and instrument placement procedures. The Multi-Pro 2000 is widely used for ultrasound-guided fine-needle and core biopsy procedures, fluid aspiration and catheterization procedures. The Multi-Pro 2000 includes inserts for 14, 16, 18, 20, 22 and 25 gauge instruments. Select kits are conveniently packaged with gel and bands. U.S. Patent No. 5,941,889.

REF. #	Product Description	Qty.
610-131	Sterile Multi-Pro 2000	24
610-133	Sterile Multi-Pro 2000 with 14 x 61 cm (5.5" x 24") CIV-Flex cover (3D)	24
610-134	Sterile Multi-Pro 2000 with 17.8 x 61 cm (7" x 24") CIV-Flex cover	24
610-137	Sterile Multi-Pro 2000 with 8 x 45cm (3.2" x 18") latex cover	24
610-138	Sterile Multi-Pro 2000 with 10 x 45cm (3.9" x 18") latex cover	24

Bulk Packaged Latex-Free CIV-Flex™ Covers



Bulk packaged in an easy-to-use box, these non-sterile CIV-Flex covers offer added protection against cross-contamination between routine scanning procedures. Non-sterile colored elastic bands are included. These covers offer a three-dimensional "box end."

REF. #	Product Description	Qty.
610-362	Non-sterile 14 x 61 cm (5.5" x 24") flat-folded CIV-Flex cover (3D)	50
610-903	Non-sterile 14 x 30cm (5.5" x 11.8") flat-folded CIV-Flex cover (3D)	50

Latex-Free CIV-Flex™ Covers



NEW Cover Sizes for Small Footprint Transducers!

Recommended for over a decade as the latex-free cover of choice, CIV-Flex covers are now available in extended lengths of 36 and 58 inches. These extended length covers offer additional sterile protection when performing routine puncture and drainage procedures. The improved telescopic-fold allows for easy application of gel, as well as easy application of the cover onto the transducer. Sterile general purpose procedure kits include transducer cover, gel packet and colored elastic bands. Select covers offer a three-dimensional "box end."

REF. #	Product Description	Qty.
610-323	Sterile 6 x 91.5cm (2.5" x 36") flat-folded CIV-Flex cover	24
610-001	Sterile 8.9 x 61cm (3.5" x 24") flat-folded CIV-Flex cover	24
610-637	Sterile 8.9 x 91.5cm (3.5" x 36") telescopically-folded CIV-Flex cover	24
610-1000	Sterile 10.2 x 147cm (4" x 58") telescopically-folded CIV-Flex cover	24
610-002	Sterile 14 x 61cm (5.5" x 24") flat-folded CIV-Flex cover (3D)	24
610-542	Sterile 14 x 91.5cm (5.5" x 36") telescopically-folded CIV-Flex cover (3D)	24
610-575	Sterile 14 x 147cm (5.5" x 58") telescopically-folded CIV-Flex cover (3D)	24
610-004	Sterile 17.8 x 61cm (7" x 24") flat-folded CIV-Flex cover	24
610-1010	Sterile 17.8 x 147cm (7" x 58") telescopically-folded CIV-Flex cover (3D)	24
610-1035	Sterile 17.8 x 147cm (7" x 58") telescopically-folded CIV-Flex cover (3D) w/ tape strips (no gel or bands)	24
610-005	Sterile 22.9 x 61cm (9" x 24") flat-folded CIV-Flex cover	24

Latex Covers



Latex covers are rolled for easy application to a variety of curved array transducers. An extended length in the larger sizes of latex offers additional sterile coverage during general puncture and drainage procedures. Latex covers are packaged with sterile gel packet and colored elastic bands.

REF. #	Product Description	Qty.
610-042	Sterile 4 x 30 (1.6" x 11.8") rolled latex cover	24
610-043	Sterile 5 x 30cm (2" x 11.8") rolled latex cover	24
610-044	Sterile 8 x 45cm (3.2" x 18") rolled latex cover	24
610-046	Sterile 10 x 45cm (3.9" x 18") rolled latex cover	24

Latex-Free Polyethylene Covers



NEW Extended Length Covers for Added Protection!

Economically-priced polyethylene covers are commonly used in sterile scanning and puncture procedures. Polyethylene covers are a safe means for reducing cross-contamination. Polyethylene covers are packaged with sterile gel packet and elastic bands. Select covers offer a three-dimensional "box end."

REF. #	Product Description	Qty.
610-013	Sterile 5 x 61cm (2" x 24") flat-folded polyethylene cover	24
610-015	Sterile 10.2 x 61cm (4" x 24") flat-folded polyethylene cover (3D)	24
610-208	Sterile 15.3 x 83.8cm (6" x 33") flat-folded polyethylene cover	24
610-020	Sterile 17.8 x 61cm (7" x 24") flat-folded polyethylene cover	24
610-022	Sterile 22.9 x 61cm (9" x 24") flat-folded polyethylene cover	24
610-989	Sterile 15.3 x 147cm (6" x 58") telescopically-folded polyethylene cover	24

Visit www.civco.com for photos, pricing and product spec sheets.

Disposable Needle Guides

- Reduces technique variability, reducing procedure time
- Contoured design for patient comfort
- Disposable design reduces risks of cross-contamination
- Most guides accept 16-18 gauge instruments
- Packaged with cover, gel and bands to increase efficiencies

Reusable Needle Guides

- Reduces technique variability, reducing procedure time
- Typically manufactured from stainless steel allowing autoclave sterilization
- Sleek design provides added patient comfort
- Most guides accept 16-18 gauge instruments

EC4-9ES • EC4-9IS • ER4-9ES/140D • 10R EV



Disposable Transrectal Needle Guide

REF. #	Product Description	Gauge Size	Qty.
657-014	Sterile transrectal needle guide	16-18	24
610-543	Sterile transrectal needle guide, 2 x 20cm (.8" x 7.9") latex cover	16-18	24

EC4-9ES • EC4-9IS • ER4-9ES/140D



Reusable Endocavity Needle Guide

REF. #	Product Description	Gauge Size	Qty.
610-666	Non-sterile endocavity needle guide with 2 x 20cm (.8" x 7.9") latex covers (10), 3.5 x 20cm (1.4" x 7.9") latex covers (10), cleaning brush (5)	16-18	1
610-1025	Non-sterile endocavity thumbscrew (Ref. # 079-400)		5

EC4-9/10ED • EC4-9/13CD



Reusable Endocavity Needle Guide

REF. #	Product Description	Gauge Size	Qty.
610-1012	Non-sterile endocavity needle guide, with cleaning brush (5)	16-18	1
610-1025	Non-sterile endocavity thumbscrew (Ref. # 079-400)		5

NEV4-9ES • NER4-9ES



NEW! Disposable Endocavity Needle Guide

REF. #	Product Description	Gauge Size	Qty.
610-958	Sterile endocavity needle guide	16-18	24
610-959	Sterile endocavity needle guide with 2.6 x 20cm (.8" x 7.9") latex cover, 3.5 x 20cm (1.4" x 7.9") latex cover	16-18	24

NEV4-9ES • NER4-9ES



NEW! Reusable Endocavity Needle Guide

REF. #	Product Description	Gauge Size	Qty.
657-024	Non-sterile endocavity needle guide with cleaning brush (5)	16-18	1
610-1025	Non-sterile endocavity thumbscrew (Ref. # 079-400)		5

3D4-9ES



Reusable Endocavity Needle Guide

REF. #	Product Description	Gauge Size	Qty.
657-025	Non-sterile 3D4-9ES endocavity needle guide with cleaning brush (5)	16-22	1
610-1025	Non-sterile endocavity thumbscrew (Ref. # 079-400)		5

Latex-Free NeoGuard™ Covers



NeoGuard is a rolled, latex-free material specially designed for ultrasound use. Recommended for patients and professionals with latex protein sensitivities, these high-quality covers are available for popular endocavity transducers. Covers include colored elastic bands. Sterile gel packets are available separately.

REF. #	Product Description	Qty.
610-843	Sterile 2.6 x 30cm (1" x 11.8") rolled NeoGuard cover	24
610-1038	Sterile 2.6 x 30cm (1" x 11.8") rolled NeoGuard cover	50
610-836	Non-sterile 2.6 x 30cm (1" x 11.8") rolled NeoGuard cover	24
610-1039	Non-sterile 2.6 x 30cm (1" x 11.8") rolled NeoGuard cover	50
610-844	Sterile 4 x 30cm (1.6" x 11.8") rolled NeoGuard cover	24
610-838	Non-sterile 4 x 30cm (1.6" x 11.8") rolled NeoGuard cover	24

Latex-Free CIV-Flex™ Covers



Especially popular for OB/GYN scanning procedures, CIV-Flex is a soft, pliable material comfortable for patients. The durable CIV-Flex material is acoustically transparent, offering distortion-free imaging. Designed specifically for ultrasound use, CIV-Flex slides easily over transducers and is a recommended alternative for latex-sensitive patients and professionals. Sterile cover kits include gel packet and colored elastic bands.

REF. #	Product Description	Qty.
610-006	Sterile 11.9 tapered to 3.8 x 61cm (4.7" tapered to 1.5" x 24") CIV-Flex cover	24
610-007	Non-sterile 11.9 tapered to 3.8 x 61cm (4.7" tapered to 1.5" x 24") CIV-Flex cover (no gel)	24

Latex PROcovers™



These high-quality, custom-sized PROcovers provide a defense against cross-contamination and protect endocavity transducers. All covers include colored elastic bands.

REF. #	Product Description	Qty.
610-075	Sterile 2 x 20cm (.8" x 7.9") rolled latex cover	24
610-039	Non-sterile 2 x 20cm (.8" x 7.9") rolled latex cover	24
610-212	Sterile 2 x 30cm (.8" x 11.8") rolled latex cover	24
610-076	Sterile 2.6 x 20cm (1" x 7.9") rolled latex cover	24
610-040	Non-sterile 2.6 x 20cm (1" x 7.9") rolled latex cover	24
610-213	Sterile 2.6 x 30cm (1" x 11.8") rolled latex cover	24
610-101	Non-sterile 2.6 x 30cm (1" x 11.8") rolled latex cover	24
610-042	Sterile 4 x 30cm (1.6" x 11.8") rolled latex cover, gel	24
610-041	Non-sterile 4 x 30cm (1.6" x 11.8") rolled latex cover	24

3.5 x 20cm Latex Covers



Conveniently packaged in individual perforated packets, the 3.5 x 20cm, a looser fitting rolled latex cover, provides protection from cross-contamination in all types of endocavity applications. Colored elastic bands are included with each cover.

REF. #	Product Description	Qty.
610-214	Sterile 3.5 x 20cm (1.4" x 7.9") rolled latex cover	50
610-011	Non-sterile 3.5 x 20cm (1.4" x 7.9") rolled latex cover	24
610-010	Non-sterile 3.5 x 20cm (1.4" x 7.9") rolled latex cover	50

Bulk Packaged Latex Covers



These cost-effective transducer covers are bulk packaged in protective light-sensitive bags for extended shelf life. Economically priced as low as 17 cents per cover, the covers provide protection against cross-contamination. Color elastic bands not included, but are sold separately.

REF. #	Product Description	Qty.
610-249	Non-sterile 3.5 x 20cm (1.4" x 7.9") rolled latex cover	100
610-262	Non-sterile 3.5 x 20cm (1.4" x 7.9") rolled latex cover	500
610-690	Non-sterile .8cm (.31") latex bands, small, green	500
610-691	Non-sterile latex bands, .8cm (.31"), small, green	10,000

Eclipse® Covers



The Eclipse transducer cover is a pre-gelled, latex-free endocavity transducer cover. The cover is made of a pliable, polyethylene material. Use of the cover reduces the concerns of cross-contamination and protects the transducer during ultrasound procedures. Each non-sterile cover is pre-gelled with Aquasonic® 100 transmission gel and is individually packaged.

REF. #	Product Description	Qty.
610-792	Non-sterile 6.3 tapered to 4.6 x 24.1 cm (2.5" tapered to 1.8" x 9.5") polyethylene cover with gel	100

Soaking System for Endocavity Transducers



This newly designed lightweight soaking cup provides a safe and efficient soaking system for transvaginal and transrectal ultrasound transducers. This product now offers a larger blue funnel for use with larger style endocavity transducers. The soaking cup kit includes a 12" x 5" high-density polyethylene cup, lid, mounting hardware and contoured funnels for smaller transducers with a width of 2.2cm (.88") and a new funnel for larger transducers with a width of 2.7cm (1.08"). The larger funnel is now colored blue for easy identification. The contoured funnel prevents transducers from hitting the bottom of the cup, saving costly damage to the lens of the transducer. The soaking cup holds up to 120 fluid ounces of disinfecting liquid and will accommodate two endocavity transducers measuring up to 22.9cm (9") in length. When not in use, the lid is provided to prevent unpleasant fumes and evaporation of disinfectant or sterilizing solutions. The cup is recommended for use with glutaraldehyde cold chemical disinfecting/sterilizing solutions, including the Cidex® and MetriCide® families of disinfecting solutions. It is recommended the cup, lid and funnels are replaced when the disinfectant expires. Replacement cups, lids and funnels are sold separately.

NEW
3D Funnel
Available!

REF. #	Product Description	Qty.
610-584	Endocavity soaking cup kit with polyethylene cup, lid, funnel and mounting hardware	1
610-585	Replacement kit with cup, lid and funnel	3
610-586	Soaking cup stabilizer, rubber feet (5), assembly hardware and stabilizer bracket	1

Intraoperative Cord Covers & Drapes



Cord covers and system drapes provide a wide sterile field during ultrasound procedures in the operating room or surgical suite. Probe clips and elastic bands are included with all cord covers. Adhesive tape is included with each system drape.

REF. #	Product Description	Qty.
610-024	Sterile 7.6 x 244cm (3" x 96") polyethylene cord cover	24
610-025	Sterile 76 x 76cm (30" x 30") polyethylene system drape	24
610-125	Non-sterile probe clips, small	100

Latex-Free Surgi™ Intraoperative Covers



The Surgi is a specially designed polyethylene transducer cover accordion-folded for easy application to transducers in a sterile environment. Providing distortion-free scanning, the Surgi extends up the length of the transducer cord, allowing up to 2.4 meters (eight feet) of sterile protection during ultrasound procedures in the operating room or surgical suite. Select covers include sterile gel packets and bands. The Surgi also is available flat-folded.

REF. #	Product Description	Qty.
610-023	Sterile 15.2 tapered to 7.6 x 244cm (6" tapered to 3" x 96") accordion-folded cover (no gel)	12
610-367	Sterile 15.2 tapered to 7.6 x 244cm (6" tapered to 3" x 96") accordion-folded cover	12
610-217	Sterile 15.2 tapered to 7.6 x 244cm (6" tapered to 3" x 96") flat-folded cover	24

Latex & Latex-Free Surgi-Tip™ Intraoperative Covers



The Surgi-Tip family of intraoperative covers offers a snug fit to small footprint intraoperative transducers. Available with latex or latex-free NeoGuard™ tips, the covers are accordion or telescopically-folded for easy application in a sterile environment. Surgi-Tip covers provide extended sterile coverage for your transducer and cord. Covers include sterile gel packet, probe clips and colored elastic bands.

REF. #	Product Description	Qty.
610-657	Sterile 15.2 x 244cm (6" x 96") poly accordion-folded cover with attached 3cm (1.2") latex tip with clips	12
610-833	Sterile 15.2 x 244cm (6" x 96") poly telescopically-folded cover with attached 3cm (1.2") latex-free NeoGuard tip with clips	12
610-776	Sterile 20.3 x 244cm (8" x 96") poly telescopically-folded cover with attached 6cm (2.4") latex tip with clips	12
610-834	Sterile 20.3 x 244cm (8" x 96") poly telescopically-folded cover with attached 6cm (2.4") latex-free NeoGuard tip with clips	12

Foam Tip Protector



These single-use cylindrical foam protectors fit over the distal tip of laparoscopic and pediatric transesophageal transducers to protect the transducer from damage during storage and transportation between patients.

REF. #	Product Description	Qty.
610-748	Non-sterile foam tip guards	24

Visit www.civco.com for photos, pricing and product spec sheets.

Latex-Free CIV-Flex™ TEE Covers & Bite Blocks



NEW! IMPROVED PRODUCT • This latex-free transesophageal echocardiography (TEE) cover features a custom tapered design to provide a tighter fit to TEE transducers. The new improved cover provides a defense against cross-contamination, while protecting the transducer. The form-fitting cover is pre loaded on the applicator for ease of application to TEE transducers. Bite blocks provide a safe passage through the patient's mouth. Select kits include bite block and gel-filled syringe.

REF. #	Product Description	Qty.
610-979	Non-sterile latex-free E-Z Guard™ pediatric bite block with strap	24
610-160	Non-sterile latex-free E-Z Guard™ adult bite blocks without strap	24
610-749	Non-sterile latex-free E-Z Guard™ adult bite blocks with strap	24
610-933	Sterile 4.6 tapered to 2.5 x 122cm (1.8" tapered to 1" x 48") CIV-Flex cover on applicator, gel packet	12
610-840	Non-sterile 4.6 tapered to 2.5 x 122cm (1.8" tapered to 1" x 48") CIV-Flex cover on applicator, gel-filled syringe, bite block	12
610-860	Non-sterile 4.6 tapered to 2.5 x 122cm (1.8" tapered to 1" x 48") CIV-Flex cover on applicator	12

TEE Fume-Free Soaking Station



This unique curved soaking tube design makes it easy to insert and remove long TEE transducers. The heavy-duty stainless steel safety cage protects against damage to the tubes and probes. The soaking station is available with two or three tubes. One tube is used for initial rinsing in clean water. Each tube holds less than one quart of disinfectant, as compared to two gallons used in conventional soaking bowls. This amounts to savings of 45.5 gallons of disinfectant a year. This soaking system meets all OSHA and JCAHO guidelines. Contact PCI Medical at 800-862-3394 for replacement filters. Use product reference number F14 for reorder. Dimensions: 48" x 14" x 15"

REF. #	Product Description	Qty.
610-850	Soaking station for TEE transducers, 110-220V (Ref. # G14-TC2)	1
610-851	Wall-mounted storage tube (Ref. # AC-SUTC)	1

General Purpose Standoffs

	REF. #	QTY.
Disposable Aquaflex gel pads (Ref.# 04-02)	610-113	6
Acoustic standoff, 1 x 10 x 15cm (.4" x 4" x 6")	610-459	1
Acoustic standoff, 2 x 10 x 15cm (.8" x 4" x 6")	610-460	1
Sterile Hydrosan couplant sheets, 10 x 15cm (4" x 6")	610-346	5

Ultrasound Gel

	REF. #	QTY.
Aquagel® 5 oz. flip-top tube	610-920	20
Aquasonic® 100, 20 gram sterile gel packets (Ref. #01-01)	610-117	48
Aquasonic® 100, .25 liter bottles (Ref. #01-08)	610-130	12
Aquasonic® 100, 5 liters, 1 dispenser bottle (Ref. #01-50)	610-116	1
Aquasonic® Clear, 5 liters econopack, 2 dispenser bottles (Ref. #03-54)	610-455	4
Scan™, 237ml bottle (Ref. # 11-08)	610-892	12
Scan™, 1 gallon pump, 2 dispenser bottles (Ref. #11-28S)	610-123	4
Ultra/Phonic, 20ml sterile gel packets (Ref. #NDC36-1000-20)	610-057	25

<i>Gel Warmers</i>	<i>REF. #</i>	<i>QTY.</i>
Gel and lotion warmer, multi-bottle-110V (Ref. # 201-786)	610-993	1
NEW • Thermasonic®, single-bottle gel warmer-120V (Ref. #82-01)	610-1020	1
NEW • Thermasonic®, single-bottle gel warmer-230V (Ref. 82-01-20)	610-1021	1
NEW • Thermasonic®, 3-bottle gel warmer-120V (Ref. #82-03)	610-1022	1
NEW • Thermasonic®, 3-bottle gel warmer-230V (Ref.82-03-20)	610-1023	1
<i>Transducer Positioning</i>	<i>REF. #</i>	<i>QTY.</i>
Non-sterile Positioning Arm, 32.4" (82.3cm)	610-200	1
Ultrasound Transducer Holder- Medium Handle	610-203	1
Rail Adapter, 1.1" x .38" (2.8 x 1cm)	610-101	1
Table Adapter, 1" (2.5cm)	610-100	1
Sterile 10" tapered to 2" x 58" (25.4 tapered to 5.1 x 147cm) telescopically-folded CIV-Flex Positioning Arm Cover	610-946	24
Sterile 10" tapered to 2" x 48" (25.4 tapered to 5.1 x 121.9cm) telescopically-folded CIV-Flex cover with Quick-Mount™ perforation	610-932	48
<i>Patient Positioning</i>	<i>REF. #</i>	<i>QTY.</i>
NEW • Non-sterile KanguruWeb Universal Base Unit (Ref. # RB2555)	610-1047	1
NEW • Non-sterile KanguruWeb disposable web (Ref. # DW1000)	610-1048	10
<i>Endocavity Soaking Station</i>	<i>REF. #</i>	<i>QTY.</i>
Soaking station for two endocavity transducers, 110V/60Hz (Ref. #G14KA)	610-720	1
Soaking station for one endocavity transducer, 110V/60Hz (Ref. #G10VP)	610-721	1
Starter kit with bracket, polyethylene cup, lid, smaller-white funnel, larger-blue funnel and mounting hardware	610-584	1
Replacement kit with cup, lid, smaller-white funnel (3) and larger-blue funnel (3)	610-585	3
Soaking cup stabilizer, rubber feet (5), assembly hardware	610-586	1
<i>Endocavity Accessories</i>	<i>REF. #</i>	<i>QTY.</i>
Non-sterile 2 x 14cm (.8" x 5.5") Endocavity Balloon	610-898	10
Non-sterile 2 x 7cm (.8" x 2.8") latex fingertip	610-038	100
Non-sterile 1.5 x 7cm (.6" x 2.8") latex fingertip	610-316	100
Sterile 2 x 20cm (.8" x 7.9") latex cover, 3.5 x 20cm (1.4" x 7.9") latex cover, gel, bands, o-ring and cleaning brushes, 22.9cm (9") (5)	610-667	20
Endocavity needle guide cleaning brush, 22.9cm (9"), 18GA	610-324	5
Endocavity needle guide cleaning brush, 22.9cm (9"), 18GA	610-951	200
Non-sterile latex-free o-rings, 1.3cm (.5"), black	612-052	40
Non-sterile latex-free o-rings, 1.6cm (.63"), black	610-751	100
Non-sterile latex-free o-rings, 1.9cm (.75"), white	610-145	100
Non-sterile latex-free o-rings, 2.3cm (.92"), white	610-752	25
Non-sterile latex bands, .8cm (.31"), small, green	610-690	500
<i>Scanning Kits</i>	<i>REF. #</i>	<i>QTY.</i>
Non-sterile scanning kits with extension tube, syringe and stopcock	610-224	5
Sterile extension tubes	610-256	10
Non-sterile 3.5 x 20cm (1.4" x 7.9") latex cover with 60cc syringe (1), extension tube (1), stopcock (1), bag of bands (1)	610-056	1
Non-sterile 3.5 x 20cm (1.4" x 7.9") latex cover with 60cc syringe (1), extension tube (1), stopcock (1), bag of bands (1)	610-058	1

<i>Cardiac Accessories</i>	<i>REF. #</i>	<i>QTY.</i>
Clear-backed electrode with snap (Ref. #2235)	610-738	50
Foam monitoring electrode with snap (Ref. #2259)	610-737	25
Soft cloth electrode with snap (Ref. #2238)	610-739	50
Resting electrode -adult (Ref. #CIV 310)	610-742	150
<i>TEE Positioning</i>	<i>REF. #</i>	<i>QTY.</i>
TEE probe holder for the OR	610-729	1
<i>TEE Soaking & Storage Units</i>	<i>REF. #</i>	<i>QTY.</i>
NEW • TEE transducer storage rack (Ref. #AC-TCRK)	610-1046	1
TEE wall support	610-724	1
Storage container	610-725	1
Disinfection container	610-723	1
Soaking station for TEE, 110-220V (Ref. #G14-TC2)	610-850	1
Wall-mounted storage tube (Ref. #AC-SUTC)	610-851	1
Non-sterile foam tip guards	610-748	24
<i>Disinfectants & Cleaners</i>	<i>REF. #</i>	<i>QTY.</i>
Cidex®, 14 day use, 1 gallon bottle (Ref. #2250)	610-640	4
Cidex®, 14 day use test strips, 15 strips per bottle (Ref. #2927)	610-803	2
Cidex® Plus, 28 day use, 1 gallon bottle (Ref. #2785)	610-641	4
Cidex® Plus, 28 day use test strips, 15 strips per bottle (Ref. #2926)	610-804	2
Cidex® OPA, 14 day use, 1 gallon bottle (Ref. #20390)	610-704	4
Cidex® OPA test strips, 15 strips per bottle (Ref. #20393)	610-714	2
DISPATCH® liquid spray, 22 oz. (Ref. #68967)	610-705	8
Enzol®, 1 gallon bottle (Ref. #2252)	610-643	4
NEW • Glute-Out neutralizer, 1 gallon bottle (Ref. #NG-1GAL)	610-1045	1
MetriCide 28, 28 day use, 1 quart bottle (Ref. #2805)	610-328	16
MetriTest strips, 50 strips per bottle (Ref. #304)	610-805	2
Sani-Cloth® PLUS Germicidal Disposable Wipe, (6" x 6.75") 160 sheets per dispenser (Ref. #Q89072)	610-943	12
Sani-Cloth® PLUS Germicidal Disposable Wipe, (8" x 14") 65 sheets per dispenser (Ref. #Q85084)	610-944	6
Sani-Cloth® HB Germicidal Disposable Wipe, 15.3 x 17.1cm (6" x 6.75"), 160 sheets per dispenser (Ref. #Q08472)	610-991	12
Sani-Cloth® HB Germicidal Disposable Wipe, 20.3 x 35.6cm (8" x 14"), 65 sheets per dispenser (Ref. #Q85484)	610-992	6
Towelettes, (7" x 6") 85 sheets per dispenser	610-128	6
TranseptiC® cleansing solution, 250ml pump spray (Ref. #09-25)	610-388	12
T-Spray™ transducer spray, 250ml pump spray (Ref. #NDC36-3200-25)	610-154	12
T-Spray™ II transducer spray, 250ml pump spray (Ref. #NDC36-2102-25)	610-703	12
<i>Skin Marking Pens</i>	<i>REF. #</i>	<i>QTY.</i>
Squeeze-Mark skin marking pen	610-746	20
<i>Printing Supplies</i>	<i>REF. #</i>	<i>QTY.</i>
Print frames	610-733	100
Chart insert for 1-4 110mm thermal prints	610-735	200
Chart insert for 1-12 thermal prints	610-732	50
File folder for 1-12 110mm thermal prints (paper)	610-178	50

Top-binding print holder, 1-8 thermal prints (plastic)	610-734	100
Side-binding print holder, 1-8 thermal prints (plastic)	610-736	25
Polaroid Type 667, b/w film, 10 exposures	610-481	2
Kowa KRP 10-002A, 6 rolls/box	610-761	1
Monitors	REF. #	QTY.
LMD-1950MD 19" LCD surgical monitor	610-719	1
LMD-1420MD 14" LCD medical office monitor	610-716	1
Optical Disks & CDs	REF. #	QTY.
SONY 3.5" MO Media, 128MB MOD, 512 bytes/sector (Ref. #EDM-128C)	610-630	5
SONY 3.5" MO Media, 230MB MOD, 512 bytes/sector (Ref. #EDM-230C)	610-631	5
SONY 3.5" MO Media, 540MB MOD, 512 bytes/sector (Ref. #EDM-540C2)	610-632	5
SONY 3.5" MO Media, 640MB MOD, 2048 bytes/sector (Ref. #EDM-640C2)	610-633	5
TDK 3.5" 230MB rewritable MOD	610-894	5
TDK 3.5" 640MB rewritable MOD	610-895	5
TDK 5.25" 5.2MB rewritable MOD	610-891	10
TDK rewritable CD, 700MB	610-874	25
Mitsubishi Print Media	REF. #	QTY.
NEW • CK-30S S-size printing paper and ribbon cartridge for use with CP-30DW, CP-30W and CP-31W printers.80 sheets/pack.	610-1041	3 packs/box
NEW • CK-30L L-size printing paper and ribbon cartridge for use with CP-30DW, CP-30W and CP-31W printers.50 sheets/pack.	610-1040	4 packs/box
KP-61S b/w super grade paper for P-40U, P-50U, P-60U, P-61U, P-63UM, P-65U, P-67UA, P-68U, P-90U, P-90W, P-91 SERIES, P-93W, P-93DW AND CP-700DU/DSU printers. 215 prints/roll.	610-064	4 rolls/box
KP-65HM b/w high density paper with blue background for P-40U, P-63UM, P-65U, P-67UA, P-68U, P-90U, P-90W, P-91SERIES, P-93W, P-93DW, CP-700UM, CP-700DU/DSU and CP-770DW printers. 215 prints/roll.	610-314	4 rolls/box
KP-70S b/w super grade paper for P-70U, P-71U, P-75U, P-78U and P-500W printers.	610-065	4 rolls/box
KP-91HG high gloss thermal paper for P-91W, P-91DW, P-91DW(UB), P-93W and P-93D printers.	610-812	4 rolls/box
CK-700 color paper for use with CP-700DU, CP-700U, CP-700DSU, CP-750U, CP-700UM, CP-710U and CP-770DW printers. For use with PK-700L and PK-700S ink ribbons. 200 prints/roll.	610-628	10 rolls/box
PK-700L large format printer ink ribbon for use with CP-700DU, CP-700U, CP-700DSU, CP-750U, CP-700UM, CP-710U and CP-770DW printers. 130 prints.	610-623	10 rolls/box
PK-700S standard format printer ink ribbon for use with CP-700DU/DSU, CP-700U, CP-750U, CP-700UM, CP-710U, and CP-770DW printers. 200 prints.	610-624	10 rolls/box
CK-800L color paper for CP-800UM & CP-800E series color printers. 100 prints, ink sheet included.	610-758	1 roll/box
CK-800S color paper for CP-800UM & CP-800E series color printers. 160 prints, ink sheet included.	610-759	1 roll/box
CK-900L standard paper, L-size ink sheet for CP-900UM, CP-900DW and CP-910U printers.	610-813	130 sheets/box
CK-900S standard paper, S-size ink sheet for CP-900UM, CP-900DW and CP-910U printers.	610-814	200 sheets/box
CK-980L standard paper, L-size ink sheet for CP-980DW printer.	610-930	160 sheets/box
CK-980S standard paper, S-size ink sheet for CP-980DW printer.	610-929	160 sheets/box
Mitsubishi Printers	REF. #	QTY.
Mitsubishi P-93W monochrome thermal printer	610-934	1
Mitsubishi P-93DW monochrome digital printer	610-935	1
NEW • Mitsubishi CP-910U color thermal printer	610-1044	1

Mitsubishi CP-900DW color dye sublimation printer	610-923	1
Mitsubishi CP-900UM printer	610-810	1
NEW • Mitsubishi CP-30W color printer	610-1042	1
NEW • Mitsubishi CP-30DW color printer	610-1043	1
Mitsubishi MD-3000 medical videocassette recorder	610-922	1
SONY Print Media	REF. #	QTY.
NEW • SONY DPW +120MD medical grade recording media, 5 disks per box	610-998	1
NEW • SONY MicroVault USB storage media	610-997	1
UPP-110HA b/w dynamic range paper for UP-850, UP-880, UP-890MD, UP-D890 printers. 110mm x 18m. 240 prints/roll.	610-389	10 rolls/box
UPP-110HD b/w high density paper for UP-850, UP-D860, UP-870MD, UP-880MD, UP-890MD, UP-895MD, UP-D895, UP-897MD printers. 110mm x 20m. 270 prints/roll. 5 rolls/box.	610-072	10 rolls/box
UPP-110HG b/w high density gloss paper for UP-895MD, UP-897MD printers. 110mm x 18m. 240 prints/roll. 5 rolls/box.	610-754	10 rolls/box
UPP-110S/5 (replaces UPP-110S/4) b/w standard paper for UP-810, UP-811, UP-850, UP-860, UP-870, UP-880, UP-890MD, UP-895MD, UP-D895, UP-897MD printers. 110mm x 20m. 210 prints/roll. 10 rolls/box.	610-071	10 rolls/box
UPP-210HD b/w high density paper for UP-910, UP-930, UP-960, UP-980 printers. 210mm x 25m. 130 prints/roll.	610-619	5 rolls/box
UPP-210SE b/w standard paper for UP-910, UP-930, UP-960, UP-980 printers. 210mm x 25m. 130 prints/roll.	610-620	5 rolls/box
UPT-210BL blue transparency film for use with UP-980 printer. 210mm x 12.5m. 42 prints/roll.	610-621	5 rolls/box
UPP-216HD b/w high density paper for UP-910, UP-930 printers. 216mm x 25m. 130 prints/roll.	610-320	5 rolls/box
UPC-21L color print pack (L-size paper) for the UP-20, UP-21MD, UP-D21, UP-D21MD, UP-D23MD printers. 4 ribbons, 1 cleaning sheet.	610-857	200 prints/box
UPC-21S color print pack (S-size paper) for the UP-20, UP-21MD, UP-D21MD, UP-D23MD printers. 3 ribbons and 1 cleaning sheet.	610-858	240 prints/box
UPC-510 color print pack for UP-D50, UP-51MD, UP-51MDU printers. Includes ink ribbon.	610-757	200 sheets/box
UPC-55 color print pack for the UP-D55MD printer. Includes 2 ink ribbons.	610-931	200 sheets/box
UPC-1010/1 color print pack for UP-1200, UP-1800MD, UP-1850MD, UP-D1500CN printers. Ink cartridge included.	610-395	100 prints/box
UPC-2010/3 color print pack for the UP-2100, UP-2300, UP-D2500, UP-2900MD printers. Ink cartridge included.	610-617	200 prints/box
UPC-3010 color printing pack for UP-3000, UP-3030 printers. 80 prints/box. Ink cartridge included.	610-247	80 prints/box
UPC-5010A/2 color printing pack for UP-5000, UP-5000W, UP-5200MD printers. Ink cartridge included.	610-555	100 prints/box
UPC-5020A b/w printing pack for UP-5000, UP-5200MD printers. Ink cartridge included.	610-753	100 prints/box
UPC-5510 color printing pack for UP-5500, UP-D5500, UP-5600MD, UP-5600MDU printers. Ink ribbon included for re-loading into existing cartridges.	610-390	200 prints/box
SONY Printers	REF. #	QTY.
SONY UP-55MD medical grade color A5 video printer	610-996	1
SONY UP-D897MD medical grade digital black and white printer	610-995	1
SONY DVO-1000 medical grade DVD recorder	610-994	1
SONY UP-897MD black & white printer	610-953	1
SONY UP-21MD medical grade color video printer	610-852	1
SONY UP-D23MD digital color printer	610-925	1
SONY UP-D55MD medical grade digital color A5 printer	610-926	1

SonoAce 8000SE

Transducer	Transducer Covers			Bracket Type	Biopsy Kits	
	<i>CIV-Flex™ Latex-Free Loose-Fitting</i>	<i>NeoGuard™ Latex-Free Form-Fitting</i>	<i>Latex Form-Fitting</i>		Starter Kits	Replacement Kits
C2-5ET-N	610-542		610-044			
C3-7ED-N	610-542		610-044			
C3-7EP-N	610-542		610-044			
EC4-9ES-N	610-006, 610-007	610-836, 610-843	610-010, 610-101, 610-213, 610-214	Disposable 16-18GA Reusable 16-18GA	610-666	610-543, 657-014
EC4-9/10ED-N	610-006, 610-007	610-836, 610-843	610-010, 610-101, 610-213, 610-214	Reusable 16-18GA	610-1012	
HL5-9ED	610-542		610-044			
L5-9EC-N	610-542		610-044	Reusable 2.5cm	612-085	610-608, 612-058, 612-062
L5-9EE-N	610-1010		610-046			
P2-4AH-N	610-637	610-833	610-043			

SonoAce 8000 LIVE

Transducer	Transducer Covers			Bracket Type	Biopsy Kits	
	<i>CIV-Flex™ Latex-Free Loose-Fitting</i>	<i>NeoGuard™ Latex-Free Form-Fitting</i>	<i>Latex Form-Fitting</i>		Starter Kits	Replacement Kits
2.0CW						
4.0CW						
3D3-5EK	610-1010		610-046			
3D4-7EK	610-542		610-044			
3D5-8EK	610-006, 610-007	610-836, 610-843	610-010, 610-101, 610-213, 610-214			
C2-4ES	610-637	610-833	610-043			
C2-5EL	610-542		610-044	Reusable 8cm	610-897	610-239
C2-5ET	610-542		610-044			
C3-7ED	610-542		610-044			
C3-7EP	610-542		610-044			
C4-9/10ED	610-637	610-833	610-043			
EC4-9/10ED	610-006, 610-007	610-836, 610-843	610-010, 610-101, 610-213, 610-214	Reusable 16-18GA	610-1012	
EC4-9ES	610-006, 610-007	610-836, 610-843	610-010, 610-101, 610-213, 610-214	Reusable 16-18GA Disposable 16-18GA	610-666	610-543, 657-014
HL5-9ED-N	610-542		610-044			
L5-9EC-N	610-542		610-044	Reusable 2.5cm	612-085	610-608
L5-9EE-N	610-1010		610-046			
L5-9ER-N	610-542		610-044			
P2-3AC-N	610-637	610-833	610-043			
P2-4AC-N	610-637	610-833	610-043			
P2-4AH-N	610-637	610-833	610-043			
P2-5AC-N	610-637	610-833	610-043			
P3-7AC-N	610-637	610-833	610-044			
VAW4-7	610-542		610-046			
VDW5-8B	610-006, 610-007	610-836, 610-843	610-010, 610-101, 610-213, 610-214			

SonoAce 8000EX

Transducer	Transducer Covers			Bracket Type	Biopsy Kits	
	CIV-Flex™ Latex-Free Loose-Fitting	NeoGuard™ Latex-Free Form-Fitting	Latex Form-Fitting		Starter Kits	Replacement Kits
2.0CW						
4.0CW						
C2-4ES	610-637	610-833	610-043			
C2-5EL-N	610-542		610-044	Reusable 8cm	610-897	610-239
C2-5ET-N	610-542		610-044			
C3-7ED-N	610-542		610-044			
C3-7EP-N	610-542		610-044			
C4-9/10ED-N	610-637	610-833	610-043			
EC4-9/10ED-N	610-006, 610-007	610-836, 610-843	610-010, 610-101, 610-213, 610-214	Reusable 16-18GA	610-1012	
EC4-9ES-N	610-006, 610-007	610-836, 610-843	610-010, 610-101, 610-213, 610-214	Reusable 16-18GA Disposable 16-18GA	610-666	610-543, 657-014
HL5-9ED-N	610-542		610-044			
L5-9EC-N	610-542		610-044	Reusable 2.5cm	612-085	610-608
L5-9EE-N	610-1010		610-046			
L5-9ER-N	610-542		610-044			
P2-3AC-N	610-637	610-833	610-043			
P2-4AC-N	610-637	610-833	610-043			
P2-4AH-N	610-637	610-833	610-043			
P2-5AC-N	610-637	610-833	610-043			
P3-7AC-N	610-637	610-833	610-044			

SonoAce X4

Transducer	Transducer Covers			Bracket Type	Biopsy Kits	
	CIV-Flex™ Latex-Free Loose-Fitting	NeoGuard™ Latex-Free Form-Fitting	Latex Form-Fitting		Starter Kits	Replacement Kits
C2-4ES-N	610-542		610-044			
C2-5ET-N	610-542		610-044			
C3-7ED-N	610-542		610-044			
C4-9/10ED-N	610-637	610-833	610-043			
EC4-9/10ED-N	610-006, 610-007	610-836, 610-843	610-010, 610-101, 610-213, 610-214	Disposable 16-18GA Reusable 16-18GA	610-1012	
EC4-9ES-N	610-006, 610-007	610-836, 610-843	610-010, 610-101, 610-213, 610-214	Reusable 16-18GA	610-666	610-543, 657-014
HL5-9ED-N	610-542		610-044			
L5-9EC	610-542		610-044	Reusable 2.5cm	612-085	610-608, 612-058, 612-062
L5-9EE-N	610-1010		610-046			

Visit www.civco.com for photos, pricing and product spec sheets.

SonoAce X8

Transducer	Transducer Covers			Bracket Type	Biopsy Kits	
	CIV-Flex™ Latex-Free Loose-Fitting	NeoGuard™ Latex-Free Form-Fitting	Latex Form-Fitting		Starter Kits	Replacement Kits
2.0CW						
4.0CW						
3D2-6ET	610-1010		610-046			
3D4-8ET	610-542		610-044			
3D4-8EK	610-1010		610-046			
3D5-9EK	610-006, 610-007	610-836, 610-843	610-010, 610-101, 610-213, 610-214			
C2-5EL	610-542		610-044	Reusable 8cm	610-897	610-239
C3-7EP	610-542		610-044			
HL5-12ED	610-542		610-044			
L5-12EC	610-542		610-044	Reusable 2.5cm		610-608, 612-058, 612-062
L5-12/50EP	610-542		610-044			
NER4-9ES	610-006, 610-007	610-836, 610-843	610-010, 610-101, 610-213, 610-214	Reusable 16-18GA Disposable 16-18GA	657-024	610-958, 610-959
NEV4-9ES	610-006, 610-007	610-836, 610-843	610-010, 610-101, 610-213, 610-214	Reusable 16-18GA Disposable 16-18GA	657-024	610-958, 610-959
P2-4AH	610-637	610-833	610-043			
P3-5AC	610-637	610-833	610-043			

AccuVix V10

Transducer	Transducer Covers			Bracket Type	Biopsy Kits	
	CIV-Flex™ Latex-Free Loose-Fitting	NeoGuard™ Latex-Free Form-Fitting	Latex Form-Fitting		Starter Kits	Replacement Kits
2.0CW						
4.0CW						
3D2-6ET	610-1010		610-046			
3D4-8ET	610-1010		610-046			
C2-5EL	610-542		610-044	Reusable 8cm	610-897	610-239
C2-6IC	610-542		610-044	Reusable 10cm	610-888	610-608
C3-7IM	610-542		610-044			
EC4-9IS	610-006, 610-007	610-836, 610-843	610-101, 610-213	Reusable 16-18GA Disposable 16-18GA	610-666	610-543, 657-014
L4-7EL	610-542		610-044	Reusable 2.5cm	613-203	610-133
L5-12/50EP	610-542		610-044			
L7-16IS	610-542		610-044	Reusable 1.5, 2.5, 4cm	678-002	610-608
L8-15IS	610-542		610-044	Reusable 1.5, 2.5, 4cm	678-002	610-608
P2-4AC	610-637	610-833	610-043			
P3-5AC	610-637	610-833	610-043			

Can't find your system? Refer to www.civco.com for additional system cross-reference charts.

AccuVix XQ

Transducer	Transducer Covers			Bracket Type	Biopsy Kits	
	CIV-Flex™ Latex-Free Loose-Fitting	NeoGuard™ Latex-Free Form-Fitting	Latex Form-Fitting		Starter Kits	Replacement Kits
2.0CW						
4.0CW						
3D3-5EK	610-1010		610-046			
3D4-7EK	610-1010		610-046			
3D5-8EK	610-006, 610-007	610-836, 610-843	610-010, 610-101, 610-213, 610-214			
C2-5EL	610-542		610-044	Reusable 8cm	610-897	610-239
C2-6IC	610-542		610-044	Reusable 10cm	610-888	610-608
C3-7IM	610-542		610-044			
EC4-9IS	610-006, 610-007	610-836, 610-843	610-010, 610-101, 610-213, 610-214	Reusable 16-18GA Disposable 16-18GA	610-666	610-543, 657-014
L5-9EE	610-1010		610-046			
L5-12IM	610-542		610-044			
L6-12IS	610-542		610-044	Reusable 1.5, 2.5, 4cm	678-002	610-608
L8-15IS	610-542		610-044	Reusable 1.5, 2.5, 4cm	678-002	610-608
P2-4AC	610-637	610-833	610-043			
P3-5AC	610-637	610-833	610-042			
P3-7AC	610-637	610-833	610-043			
VAW3-5	610-542		610-044			
VAW4-7	610-542		610-044			
VDW5-8B	610-006, 610-007	610-836, 610-843	610-010, 610-101, 610-213, 610-214			
VNW6-12	610-1010		610-046	Reusable 2cm	742-327	610-608

SonoAce X6

Transducer	Transducer Covers			Bracket Type	Biopsy Kits	
	CIV-Flex™ Latex-Free Loose-Fitting	NeoGuard™ Latex-Free Form-Fitting	Latex Form-Fitting		Starter Kits	Replacement Kits
2.0CW						
4.0CW						
C3-7EP	610-542		610-044			
C4-9/10ED-N	610-637	610-833	610-043			
HL5-12ED	610-542		610-044			
L5-12/50EP	610-542		610-044			
NER4-9ES	610-006, 610-007	610-836, 610-843	610-010, 610-101, 610-213, 610-214	Reusable 16-18GA Disposable 16-18GA	657-024	610-958, 610-959
NEV4-9ES	610-006, 610-007	610-836, 610-843	610-010, 610-101, 610-213, 610-214	Reusable 16-18GA Disposable 16-18GA	657-024	610-958, 610-959
P2-4AH	610-637	610-833	610-043			

Visit www.civco.com for photos, pricing and product spec sheets.

AccuVix V20

Transducer	Transducer Covers			Bracket Type	Biopsy Kits	
	CIV-Flex™ Latex-Free Loose-Fitting	NeoGuard™ Latex-Free Form-Fitting	Latex Form-Fitting		Starter Kits	Replacement Kits
2.0CW						
4.0CW						
3D2-6ET	610-1010		610-046			
3D4-8ET	610-1010		610-046			
3D4-9ES	610-006, 610-007	610-836, 610-843	610-010, 610-101, 610-213, 610-214			
3D5-9EK	610-006, 610-007	610-836, 610-843	610-010, 610-101, 610-213, 610-214			
C2-5EL	610-542		610-044	Reusable 8cm	610-897	610-239
C2-6IC	610-542		610-044	Reusable 10cm	610-888	610-608
C3-7IM	610-542		610-044			
EC4-9IS	610-006, 610-007	610-836, 610-843	610-101, 610-213	Reusable 16-18GA Disposable 16-18GA	610-666	610-543, 657-014
L4-7EL	610-542		610-044	Reusable 2.5cm	613-203	610-133
L5-12/50EP	610-542		610-044			
L5-13IS	610-542		610-044	Reusable 1.5, 2.5, 4cm	678-002	610-608
L7-16IS	610-542		610-044	Reusable 1.5, 2.5, 4cm	678-002	610-608
P2-4AC	610-637	610-833	610-043			
P3-5AC	610-637	610-833	610-043			



Worldwide Headquarters: 102 First Street South, Kalona, Iowa 52247-9589 USA
 319/656-4447, 800/445-6741 • Fax: 319/656-4451, 877/329-2482 • www.civco.com

COPYRIGHT © 2008. ALL RIGHTS RESERVED. CIVCO IS A REGISTERED TRADEMARK OF CIVCO MEDICAL SOLUTIONS. CIV-FLEX, HYDROSCAN, MAGGI, NEOGUARD, PROCOVERS, SURGI, SURGI-TIP AND ULTRA-PRO II ARE ALL TRADEMARKS OF CIVCO. ALL OTHER TRADEMARKS ARE THE PROPERTY OF THEIR RESPECTIVE OWNERS. ALL PRODUCTS MAY NOT BE LICENSED IN ACCORDANCE WITH CANADIAN LAW. CIVCO PRODUCT NAMES MAY BE COVERED BY THE FOLLOWING U.S. PATENT: 5,941,889. OTHER PATENTS PENDING. PRINTED IN USA 2008R-2377